


Specifications

REF Measurement	
Measuring range	Spherical refractive power : -30D to +25D (0.12D/0.25D steps)* Cylindrical refractive power: OD to ±12D (0.12D/0.25D steps)* Direction of astigmatic axis:0° to 180° (1°/5° steps)
Minimum measurable pupil diameter	Φ2.0mm
PD measurement range	20 to 85mm (0.5mm step)
Target fixation	Auto fog system
KRT Measurement	
Measuring range	Cornea curvature radius: 5.00mm to 13.00mm (0.01mm step) Corneal refractive power: 67.50D to 25.96D (0.12D / 0.25D steps) (where, corneal refractive power =1.3375D) Corneal astigmatic power: OD to ±12D (0.12D / 0.25D steps) Direction of corneal astigmatic axis: 0° to 180° (1°/5° steps)
Ocular Pressure Measurement	
Measuring range	1 to 60mmHg (1mmHg step)
Cornea Thickness Measurement	
Measuring range	0.400mm to 0.750mm (0.001mm step)
Others	
Chinrest travel distance	up/down: 67mm
Other Specifications	
Dimensions	293-396mm(W) x 505-601mm(D) x 470-682mm(H)
Weight	22.0kg
Power Supply	100-240V AC, 50-60Hz, 100VA

* -30D ≦ spherical refractive power + cylindrical refractive power or spherical refractive power + cylindrical refractive power ≦ +25D

Flexible data management options

measurement data can also be sent directly to your EMR system or your automated refraction head for improved efficiency.

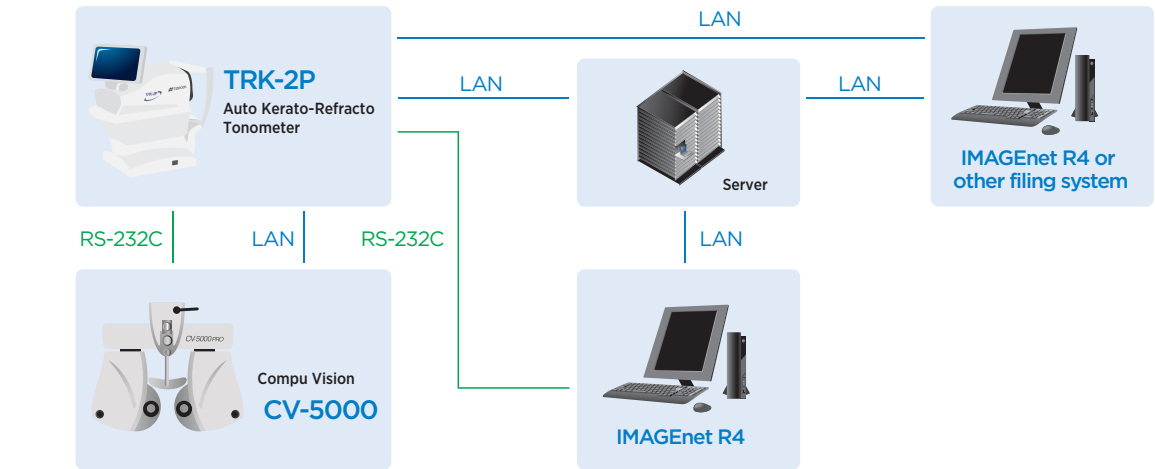

* Not available for sale in the United States.
* Not available in all countries, please check with your distributor for availability in your country
*Subject to change in design and/or specifications without advanced notice.

IMPORTANT In order to obtain the best results with this instrument, please be sure to review all user instructions prior to operation.



TOPCON CORPORATION



TOPCON CORPORATION



TOPCON CORPORATION

TOPCON CORPORATION
75-1 Hasunuma-cho, Itabashi-ku, Tokyo 174-8580, Japan. Phone:3-3558-2523/2522 Fax:3-3965-6898 www.topcon.co.jp

TOPCON MEDICAL SYSTEMS, INC.
111 Bauer Drive, Oakland, NJ 07436, U.S.A.
Phone:+1-201-599-5100 Fax:+1-201-599-5250 www.topconmedical.com

TOPCON MEDICAL LASER SYSTEMS, INC.
3130 Coronado Drive Santa Clara, California 95054 USA
Phone:+1-408-235-8200 Phone:(USA only)+1-888-760-8657
Fax:+1-408-235-8259 Email: tmlsinfo@topcon.com

TOPCON CANADA INC.
110 Provencher Avenue, Boisbriand, QC J7G 1N1 CANADA
Phone:+1-450-430-7771 Fax:+1-450-430-6457 www.topcon.ca

TOPCON EUROPE MEDICAL B.V.
Essebaan 111, 2908 LJ Capelle a/d IJssel, P.O.Box145,
2900 AD Capelle a/d IJssel, THE NETHERLANDS
Phone:+31-(0)10-4585077 Fax:+31-(0)10-4585045
E-mail: medical@topcon.nl; www.topcon.eu

ITALY OFFICE
Viale dell'Industria 80, 20037 Paderno Dugnano, (Milano), ITALY
Phone:+39-02-9186671 Fax:+39-02-91861091 E-mail: topconitaly@tiscali.it; www.topcon.it

DANMARK OFFICE
Praestemarksvej 25, 4000 Roskilde, DANMARK
Phone:+45-46-327500 Fax:+45-46-327555
E-mail: topcon@topcondanmark.dk www.topcondanmark.dk

IRELAND OFFICE
Unit 278, Blanchardstown, Corporate Park 2 Ballycoolin Dublin 15, IRELAND
Phone:+353-18975900 Fax:+353-18293915 E-mail: medical@topcon.ie; www.topcon.ie

TOPCON S.A.R.L.
BAT A1 3 route de la révolte 93206 SAINT DENIS CEDEX
Tel: +33 1 49 21 22 23 Fax: +33 1 49 21 23 24 E-mail:topcon@topcon.fr; www.topcon.fr

TOPCON DEUTSCHLAND G.m.b.H.
Hanns-Martin-Schleyer Strasse 41, D-47877 Willich, GERMANY
Phone:+49-(0)2154-8850 Fax:+49-(0)2154-885177 E-mail:med@topcon.de; www.topcon.de

TOPCON SCANDINAVIA A.B.
Neogatan 2, P.O.Box 25, 43151 Mölndal, SWEDEN
Phone:+46-(0)31-7109200 Fax:+46-(0)31-7109249 E-mail:medical@topcon.se; www.topcon.se

TOPCON ESPAÑA S.A.
HEAD OFFICE
Frederic Mompou 4 Esc. A Bajos 3, 08960 Sant Just Desvern Barcelona, SPAIN
Phone:+34-93-4734057 Fax:+34-93-4733932 E-mail: medica@topcon.es; www.topcon.es

TOPCON (GREAT BRITAIN) LTD.
Topcon House, Kennet Side, Bone Lane, Newbury, Berkshire RG14 5PX United Kingdom
Phone:+44-(0)1635-551120 Fax:+44-(0)1635-551170 E-mail:medical@topcon.co.uk; www.topcon.co.uk

TOPCON POLSKA Sp. z o.o.
ul. Warszawska 23, 42-470 Siewierz, POLAND
Phone:+48-(0)32-6705045 Fax:+48-(0)32-6713405 www.topcon-polska.pl

TOPCON SINGAPORE MEDICAL PTE. LTD.
1 JALAN KILANG TIMOR #09-01 PACIFIC TECH CENTRE SINGAPORE 159303
Phone:+65-68720606 Fax:+65-67736150 E-mail:medical_sales@topcon.com.sg www.topcon.com.sg

TOPCON INSTRUMENTS (MALAYSIA) SDN.BHD.
No. D1, (Ground Floor), Jalan Exocella 2, Off Jalan Ampang Putra,
Taman Ampang Hilir, 55100 Kuala Lumpur, MALAYSIA
Phone:+60-(0)3-42709866 Fax:+60-(0)3-42709766

TOPCON INSTRUMENTS (THAILAND) CO.,LTD.
77/162 Sinnsathorn Tower, 37th Floor, Krungthoburi Rd., Klongtongsoi,
Klongsarn, Bangkok 10600, THAILAND
Phone:+66(0)2-440-1152-7 Fax:+66-(0)2-440-1158

TOPCON CORPORATION BEIJING OFFICE
Block No.9, Kangding Street Beijing Economic-Technological Development Area,
Beijing, 100176, CHINA
Phone:+86-(0)10-6780-2799 Fax:+86-(0)10-6780-2790

TOPCON CORPORATION SHANGHAI OFFICE
14L Huamin Empire Plaza, No.726, Yan-an Xi Road,
Shanghai, 200050, CHINA
Phone:+86-(0)21-5238-7722 Fax:+86-(0)21-5237-0761

TOPCON CORPORATION BEIRUT OFFICE
P.O.Box 10-1002 Antelias Beirut, LEBANON
Phone:+961-4-523525/523526 Fax:+961-4-521119

TOPCON CORPORATION DUBAI OFFICE
P.O.Box 293705, Dubai Airport Free Zone L.L.U. J-12, Dubai, U.A.E
Phone:+971-4-299-5900 Fax:+971-4-299-5901

PRINTED WITH SOY INK
© 2013 TOPCON CORPORATION All Rights Reserved
Printed in Japan 2013 09-60 NP E101-1 GE

Auto Kerato-Refracto Tonometer
TRK-2P



TOPCON
CONNECTING VISIONS

4in1

The TRK-2P offers four times the diagnostic capabilities and is the ideal advanced patient pre-test station.

The new TRK-2P features a fully automated alignment and measurement system, combining an Auto-Refractometer, Keratometer, Non-Contact Tonometer and Pachometer in one instrument. The compact size and rotating color touch screen control panel offers total flexibility for the operator and also instrument location. The new TRK-2P improves patient flow and practice efficiency.




Features

- » Refractometer, Keratometer, Non-Contact Tonometer, and Pachometer in one single instrument
- » Fully automated measurement and alignment
- » Rotating touch screen control panel
- » Compact and modern design


Save space and patient test time

With 4 instrument functions in 1, the compact, ergonomically designed TRK-2P will free up practice space and save valuable patient test time. It can easily fit into the smallest of spaces and will help a busy practice with patient flow. Your patients will appreciate not having to be moved from instrument to another.

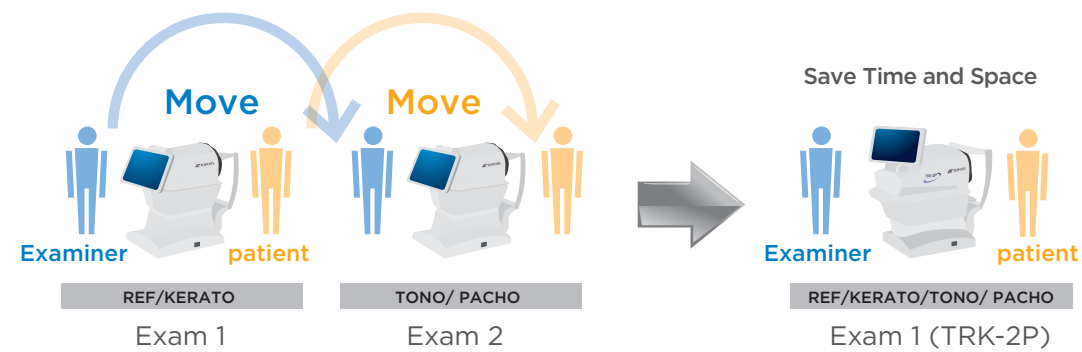

Flexible and space saving layout

The rotating touch screen control panel allows the TRK-2P to be positioned against the wall; in the corner or in the classic forward facing way, saving valuable space in the exam room. The operator also has flexibility to operate the TRK-2P from any position allowing them to observe the patient more efficiently.

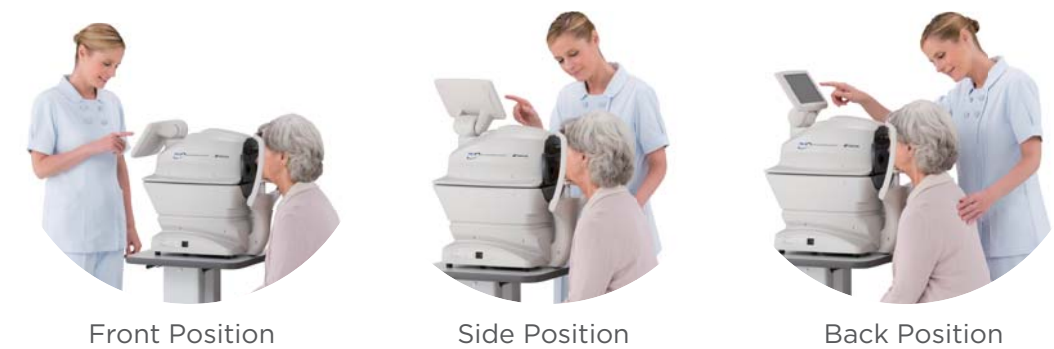


Easy-to-Use 8.5 inch color touch panel display

The large color touch panel, clearly displays the patient's eye, operating icons and all measurement results, making it very easy to use and to understand exactly how the TRK-2P is performing during the patient measurement process.

Fully automated measurement


The TRK-2P features a fully automated alignment and measurement system. With the touch of your fingertip on the color display, the TRK-2P will measure both eyes accurately and efficiently.

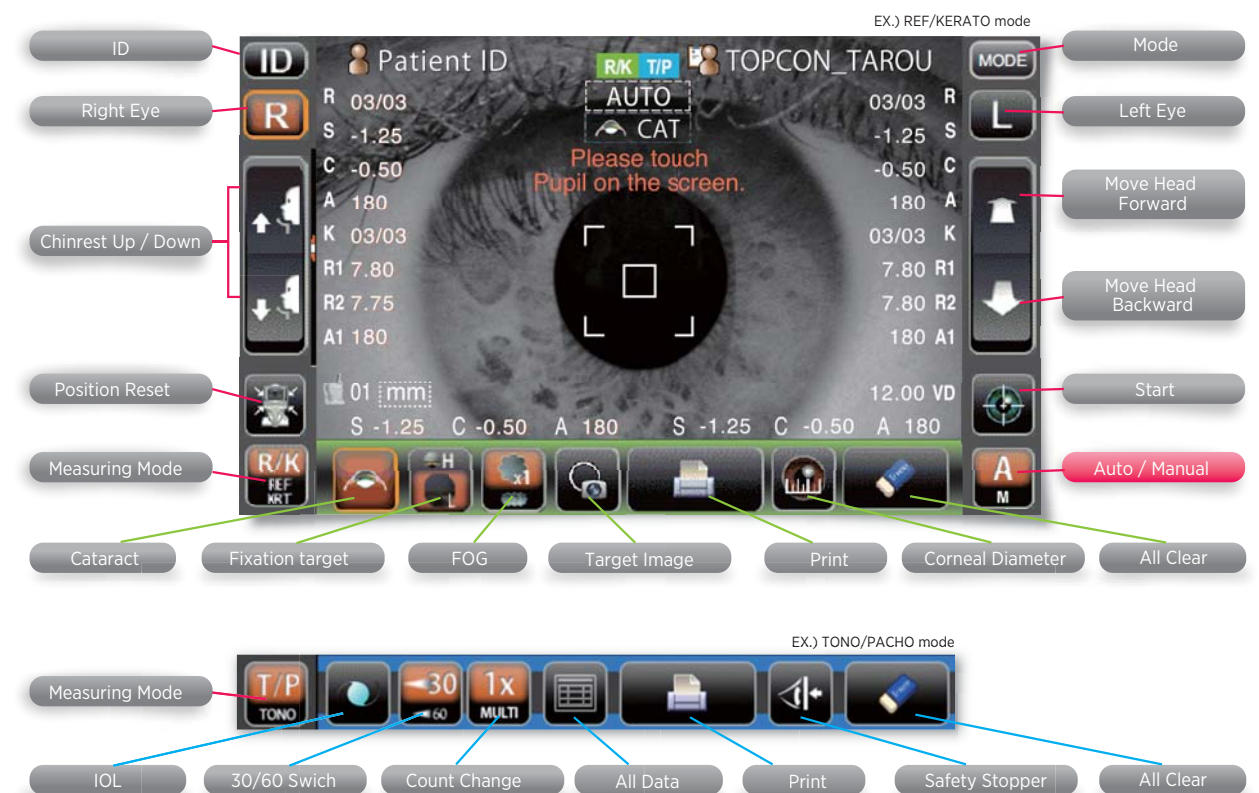

Instant trigger




For patients unable to fixate on the same point for a period of time, the TRK-2P has an instant trigger button. The operator is in full control throughout the measurement process and the trigger button can be pressed at any point, even once the auto alignment has started.

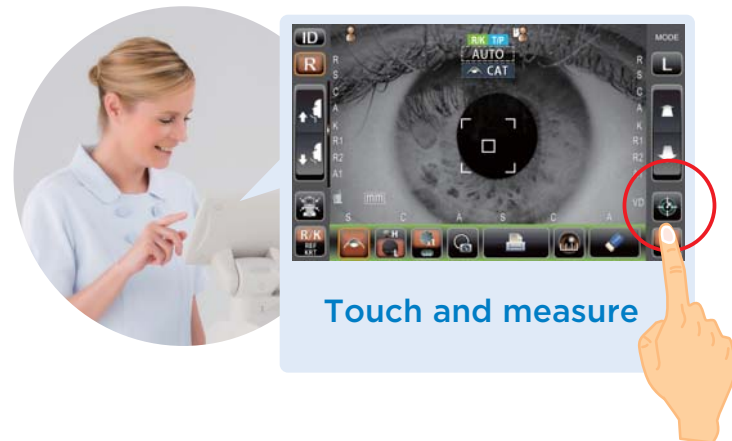

CCT and reliable (IOP) measurements




The TRK-2P provides accurate and reliable intraocular pressure measurements. Studies show that the corneal thickness directly correlates to IOP measurements. The built-in pachometer function can be used to adjust the IOP value depending if the cornea is thicker or thinner than average. The integrated formula for IOP adjustment can be customised by the eyecare professional based on the latest clinical research.

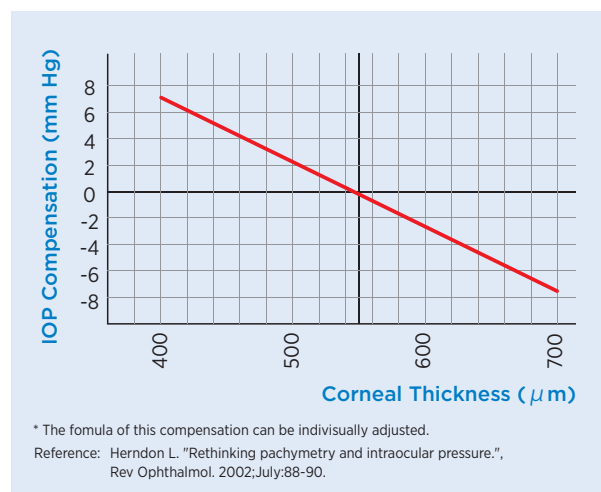

Additional Measurement Mode



If a measurement error occurs, the Additional Measurement Mode allows the TRK-2P to automatically take additional IOP measurements until all measurements are completed and error free, saving patient re-test time.

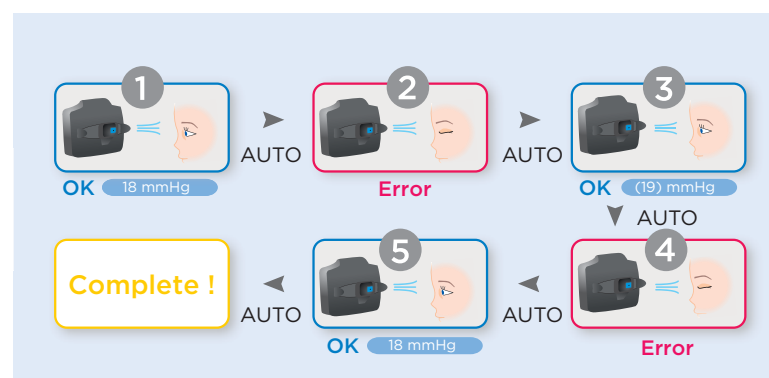

*Additional measurement mode is available in setting menu.
*Preset count can be also selected in the menu. In the example, the image shows the case preset count is 3.

Comfortable manual mode



If required, the TRK-2P can be used in manual mode. To do this, tap on the pupil image on the screen and the TRK-2P will then follow your finger position.

Accurate and reliable refractive measurements




With Topcon's exclusive Rotary Prism Technology, the TRK-2P provides unmatched accuracy and reliability. The innovative optical design incorporated into the TRK-2P allows for accurate and reliable measurements with a pupil as small as $\phi 2\text{mm}$.

Cataract mode



TRK-2P has a built-in function that allows for measurement of patients who have cataract or an opacity, this function works automatically or manually.

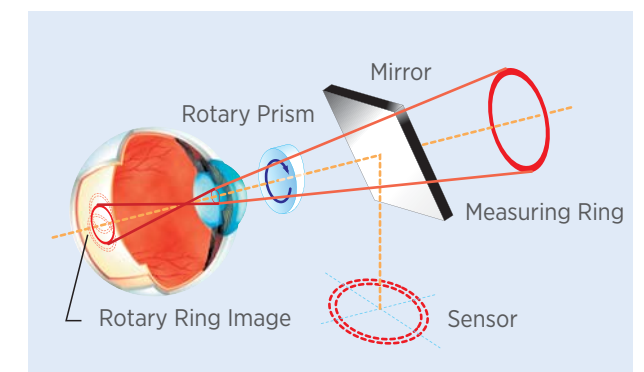

Softer air puff



Improved software controls the optimal amount of air required for every IOP measurement for each individual eye. A softer air puff improves patient comfort resulting in an improved measurement experience.

Built-in autocut printer



The TRK-2P is equipped with a built-in thermal printer that can be easily loaded with a roll of paper. The TRK-2P cuts the paper automatically at the end of printing.

